



# Medical Ozone Generator AOT-MD-500

Idea for Medical Ozone Therapy

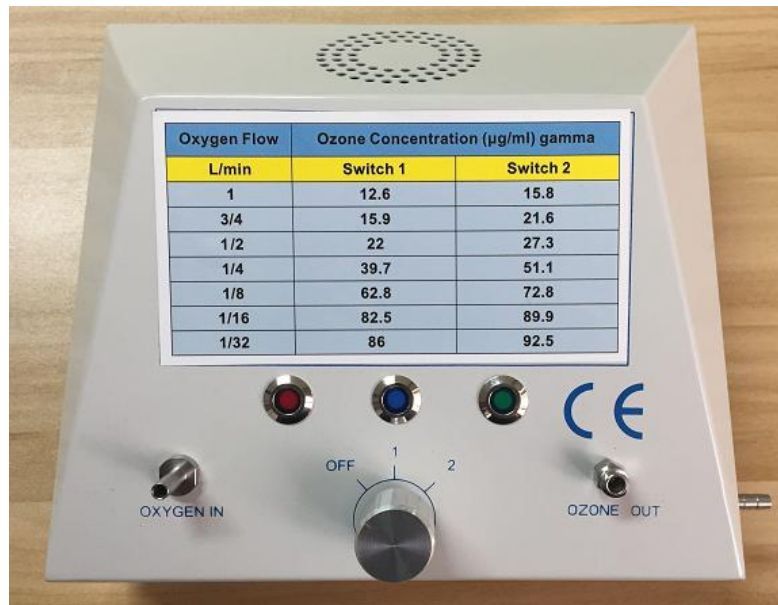

# Product introduce

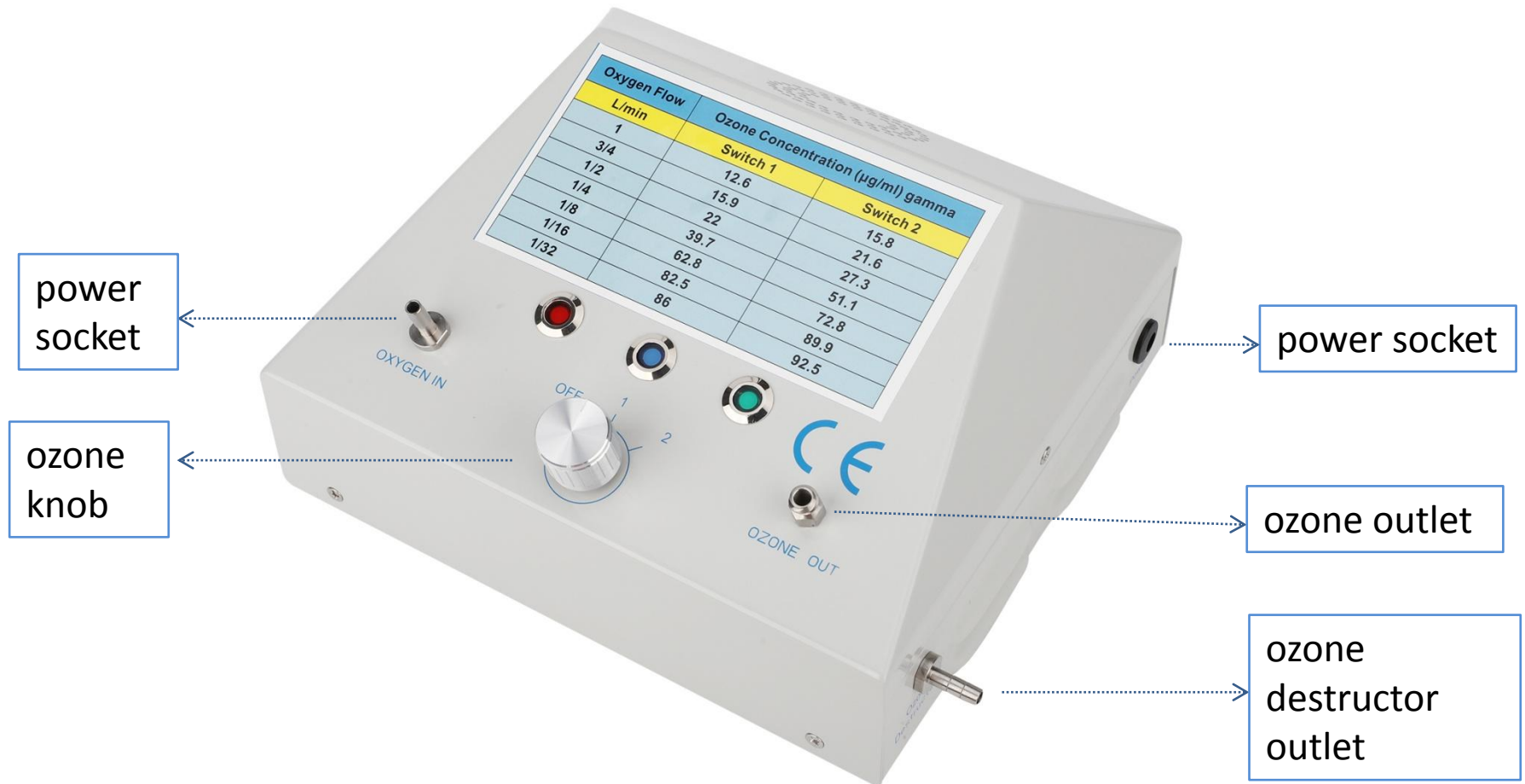

## 2 levels Ozone Concentration


| Oxygen Flow | Ozone Concentration (µg/ml) gamma |          |
|-------------|-----------------------------------|----------|
| L/min       | Switch 1                          | Switch 2 |
| 1           | 12.6                              | 15.8     |
| 3/4         | 15.9                              | 21.6     |
| 1/2         | 22                                | 27.3     |
| 1/4         | 39.7                              | 51.1     |
| 1/8         | 62.8                              | 72.8     |
| 1/16        | 82.5                              | 89.9     |
| 1/32        | 86                                | 92.5     |

# Specification



|                     |                                  |
|---------------------|----------------------------------|
| Mode No:            | AOT-MD-500                       |
| Ozone Concentration | 12.6 - 92.5 ug/ml                |
| Oxygen Flow         | 1/32L/min ----1L/min             |
| working method      | Titanium electrode & Quartz tube |
| Voltage             | AC110-220V                       |
| Power               | 15W                              |
| Product Dimension   | 206x*180* 105 mm                 |
| Net Weight          | 1.35kg                           |
| Housing Material    | Metal with powder painting       |
| Package             | 32*26.5*14 cm                    |
| Gross Weight        | 4kg                              |
| Warranty            | 1 year                           |

## key component- ozone module we use

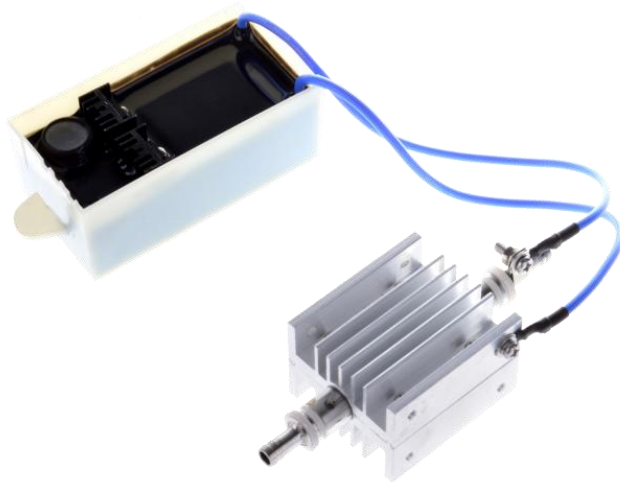

**Ozone Output** : 500mg/H @2.5 L/min,  
RH 20% ,25 °C dry air

**material** : corona discharge with Titanium  
electrode and quartz tube

**Power** : 12VDC, 10W

**Size** : 100 x 41 x 33 mm (Electronic )  
100.6 x 40 x 30 mm (Generator)

**Weight** : 0.3 kg

### Quality Guarantee:

---

**Technology** : Corona discharge

**lifespan** : 5 years

**EMC** : PASS

## Accessories Included

|                  |                                                                                      |          |
|------------------|--------------------------------------------------------------------------------------|----------|
| Luer Lock        |    | 2 pieces |
| 3-valve lock     |     | 1 piece  |
| Silicon house    |    | 2 meters |
| Diffuser stone   |    | 2 pieces |
| Adaptor 100-240V |  | 1 piece  |

## How to work ?

To run the machine, please connect with pure oxygen tank and oxygen regulator

Oxygen regulator

Oxygen Tank



you may need :

# Oxygen Regulator

AOT-OR-870

for medical type oxygen tank



AOT-OR-540

for industrial type oxygen tank



# Application



Medical ozone water



Ozone injection



Dental treatment



Major autohaemotherapy



Ear insufflation



Ozone oil



Skin therapy




Saline intravenous injection



Rectal insufflation



Pain therapy

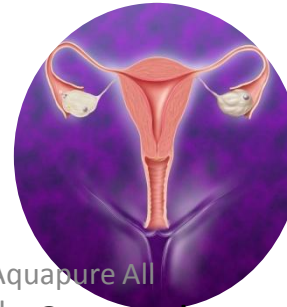

Gynecology



Ozone sauna

# Ozone Therapy Benefit

## Therapeutic Effects of Ozone therapy

1. Strong oxidation
2. Promote blood circulation
3. Accelerate hemoglobin release oxygen into tissues and improve oxygen metabolism.
4. Immunostimulation and balancing effect
5. Broad-spectrum antibacterial
6. Free radicals eradicating
7. Eliminating tumor cells

AQUAPURE is from Stuttgart Germany, the famous city with high technologies. We have famous brands of Mercedes-Benz, Porsche, Bosch, Stihl, Kaercher, etc here. It is called the 'birthplace of the automobile'.

AQUAPURE was set up in 2002. Since then we focus on health products such as plasma generator, ozone generator, ion generator technologies.

To combine the high technology from Germany and Italy with low labour cost from China, we set up AQUAURE ASIA LIMITED and AQUAPURE (Shenzhen) Ozone Technology CO.LTD for production and shipment managements in 2009 and 2015.



**Kenneth Smith**



[www.aquapureozone.com](http://www.aquapureozone.com)

Mainland China

## **AQUAPURE ASIA LTD**

R1102, University Town Business Park,  
Lishan Road, Taoyuan, Nanshan, 518055  
Shenzhen, China

**Tell :** 0086-755- 26605963

Erika.wang@aquapureasia.com

Germany

## **AQUAPURE EUROPE GmbH**

Hans-Thoma-Weg 3, 71032 Boeblingen, Germany

Tell: 0049 (0) 7031 22 96 88

Mobil: 0049 (0) 1715 21076

kenneth@aquapuregroup.eu